

SP-1P

Specular microscope



Functionality & aesthetics combined



The next generation in specular microscopy

The Topcon specular microscope SP-1P introduces a modern ergonomic design and innovative features that simplifies its use and increases its efficiency. A large 10.4 inch rotatable touchscreen eliminates the need for a control lever, and can be turned to various angles for better interaction with the patient.

By simply tapping on the center of the patients' pupil displayed on the monitor, the SP-1P automatically centers, focuses, acquires and analyzes the endothelial cell image. The entire operation takes a few seconds and it is fast and smooth, requiring minimum training.

The SP-1P offers space saving features and has a very easy to use interface.

Features

Compact and stylish design

- 10.4 inch rotatable touchscreen

Quick automatic measurement and analysis

- Instant acquisition of the analysis results
- Intuitive operation

Panorama photography function

- Substantial size increase of the analyzed area

Easy report

- Built-in printer
- Export to storage

Two photography modes

- Sequence mode
- Free style mode

Easy to read screen and comprehensive analysis software

- Frequently referred values are shown on top
- A pleomorphic / polymorphic histogram can be shown with color



Features



Compact and stylish design


Topcon specular microscope SP-1P offers comfort for patients, ease of use, and increased operation speed. The SP-1P can be placed against the wall, because of the rotatable touchscreen.



Quick automatic measurement and analysis


The SP-1P photographs and displays analyzed data by simply tapping the screen. All functions, alignment, photography and analysis are quickly performed automatically.

Conventional analysis area

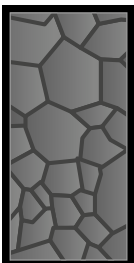

Eg. N 40

Panorama

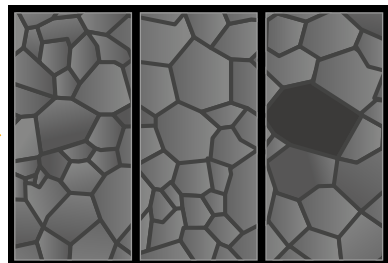

Eg. N 100


Panorama photography function

The 'panorama' function takes three images in three different areas: central, nasal and temporal and automatically combines them creating a larger area for the observation and analysis of endothelial cells*. Compared to conventional analysis, this larger area of cells allows for a faster and more complete evaluation of the patients' endothelium condition. The 'panorama' function is especially useful when the cells of the patient are decreased and the available number of cells is limited. The increased analyzed cells improve on the reliability of the examination result.

* The analysis of combined images is available after photographing a minimum of 2 areas (center and either adjacent nasal / temporal).

Features

Single

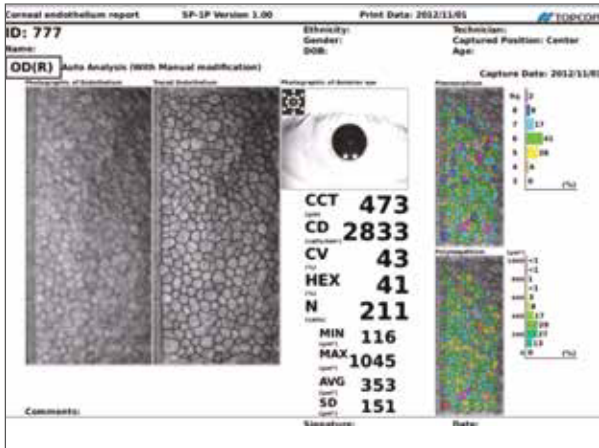


Easy report


Sample report

After photography the report will be transferred to the selected storage location.


Panorama

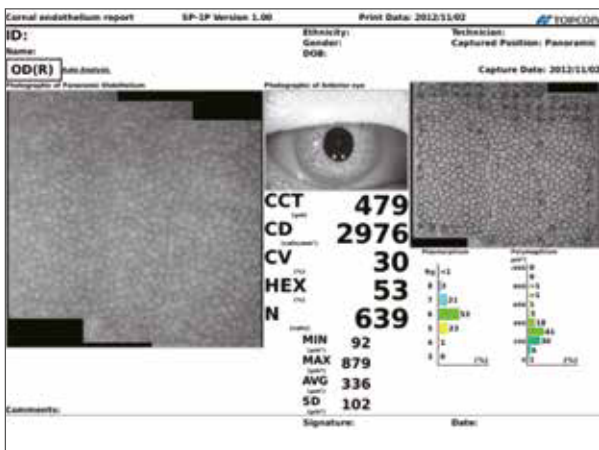

Single

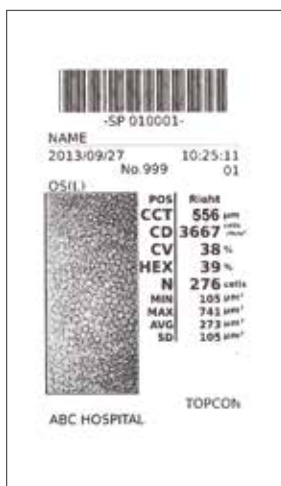

Panorama

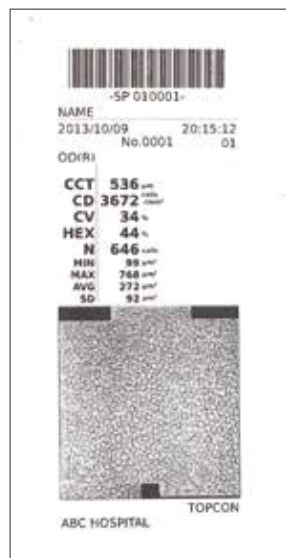

Sample printout

Direct printout from the device.

Features

Two photography modes

The SP-IP offers two photography modes. The examiner can select between sequence mode and free style mode.

Sequence mode

The examiner can enter patient information prior to photographing and then perform the photograph

and analysis of the different areas of the cornea. This ensures that all necessary information on the patient is saved.




All operation can be done by the touchscreen, simply by following the instructions.



Patient registration




Select the eye and location to examine

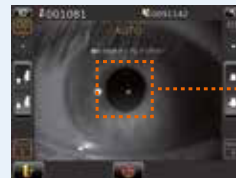

Central Photo

Tap the center of the pupil

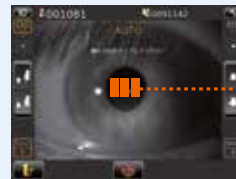

Panorama

Center/Nasal/Temporal positions can be photographed and integrated



Peripheral Photo

Photograph peripheral area(s)



Power on




Start examination


Free style mode

Immediately after turning on the device, the examiner will be able to take pictures. Even during the examination, the examiner can change or add the area to be photographed which allows a more flexible and quick analysis.

Easy to read screen

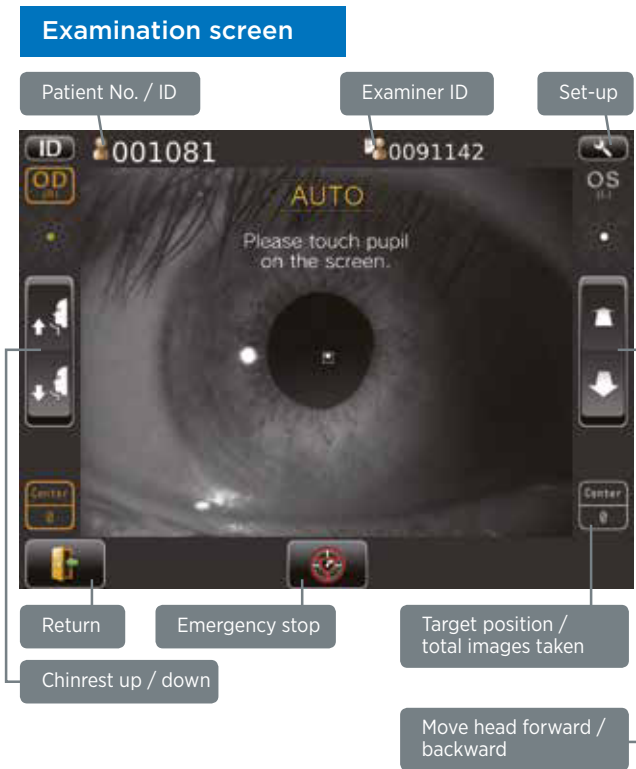

Display analyzed result

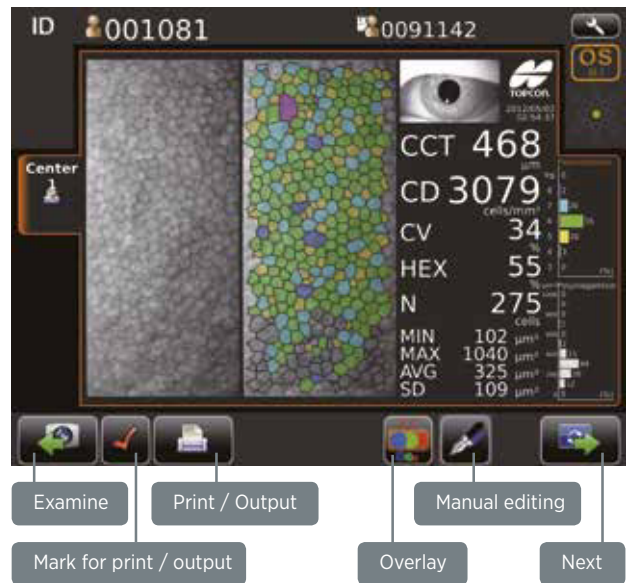


Comprehensive analysis software


All relevant values such as corneal thickness (CCT), cell density (CD), coefficient of variation (CV) and percentage of cell hexagonality (HEX) are displayed

in large font size on top of the measurement screen. By tapping the overlay button at the bottom, the analysis graph can be changed to a pleomorphic or polymegathic histogram.

Single

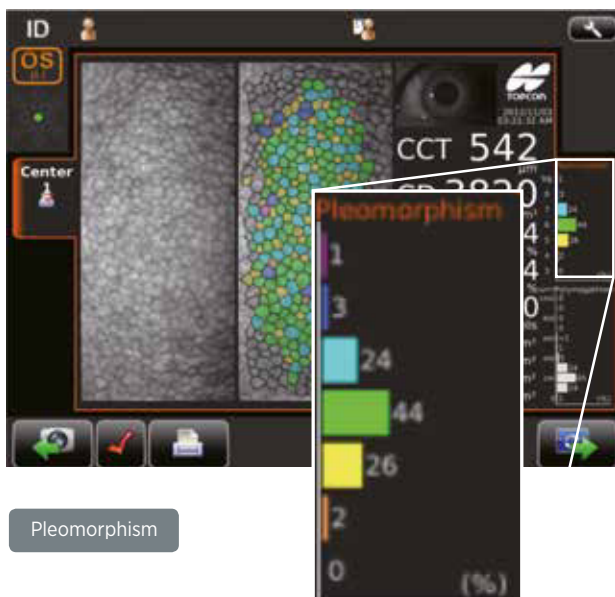

Panorama

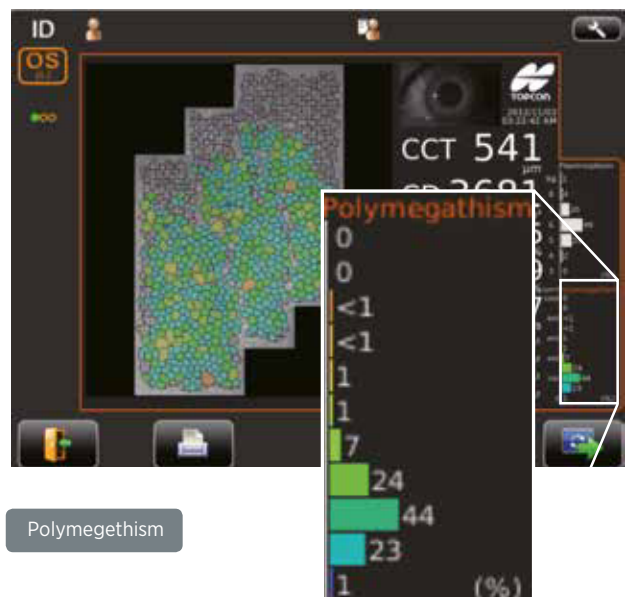


Specifications

| | |
|----------------------------------------|---------------------------------------------------------|
| Corneal endothelium photography | |
| Photography magnification | 254x (on the control panel) |
| Photography range | 0.25 mm x 0.55 mm |
| Revolving power | More than 125 line/mm |
| Fixation target | Central and (semi) peripheral |
| Corneal thickness measurement | |
| Measurement range | 0.400 - 0.750 mm Display unit: 0.001 mm step display |
| Other specifications | |
| Dimensions | 286 - 468 mm (W) x 445 - 592 mm (D) x 486 - 681 mm (H) |
| Weight | 17 kg |
| Power supply | Source voltage |
| | Power input |
| | 100-240V AC, 50-60Hz |
| | 70-120VA |

* Corneal thickness is measured only when taking a photograph using the central fixation target.

System chart

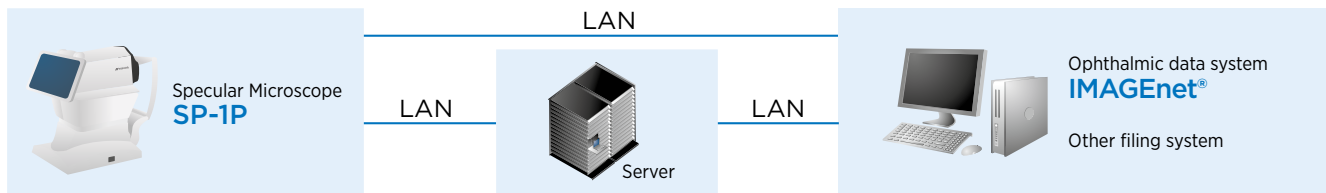

IMPORTANT

Subject to change in design and/or specifications without advanced notice.
In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.
Medical device Class Im. Manufacturer: Topcon Corporation.



TOPCON CORPORATION

TOPCON CORPORATION

Topcon Europe Medical B.V.
Essebaan 11; 2908 LJ Capelle a/d IJssel; P.O. Box 145;
2900 AC Capelle a/d IJssel; The Netherlands
Phone: +31-(0)10-4585077; Fax: +31-(0)10-4585045
E-mail: medical@topcon.eu; www.topcon-medical.eu

Topcon Danmark
Præstemarksvej 25; 4000 Roskilde, Denmark
Phone: +45-46-327500; Fax: +45-46-327555
E-mail: info@topcon.dk
www.topcon.dk

Topcon Scandinavia A.B.
Neogatan 2; P.O. Box 25; 43151 Mölndal, Sweden
Phone: +46-(0)31-7109200; Fax: +46-(0)31-7109249
E-mail: medical@topcon.se; www.topcon.se

Topcon España S.A.
HEAD OFFICE; Frederic Mompou, 4;
08960 Sant Just Desvern; Barcelona, Spain
Phone: +34-93-4734057; Fax: +34-93-4733932
E-mail: medica@topcon.es; www.topcon.es

Topcon Italy
Viale dell' Industria 60;
20037 Paderno Dugnano, (MI) Italy
Phone: +39-02-9186671; Fax: +39-02-91081091
E-mail: info@topcon.it; www.topcon.it

Topcon France Medical S.A.S.
BAT A1; 3 Route de la Révolte, 93206 Saint Denis Cedex
Phone: +33-(0)1-49212323; Fax: +33-(0)1-49212324
E-mail: topconfrance@topcon.com;
www.topcon-medical.fr

Topcon Deutschland GmbH
Hanns-Martin-Schleyer Strasse 41;
D-47877 Willich, Germany
Phone: (+49) 2154-885-0; Fax: (+49) 2154-885-177
E-mail: info@topcon-medical.de; www.topcon-medical.de

Topcon Polska Sp. z o.o.
ul. Warszawska 23; 42-470 Siewierz; Poland
Phone: +48-(0)32-670-50-45; Fax: +48-(0)32-671-34-05
www.topcon-polska.pl

Topcon Great Britain Medical Ltd.
Topcon House; Kennet Side; Bone Lane; Newbury
Berkshire RG14 5PX; United Kingdom
Phone: +44-(0)1635-551120; Fax: +44-(0)1635-551170
E-mail: medical@topcon.co.uk; www.topcon.co.uk

Topcon Ireland
Unit 276, Blanchardstown; Corporate Park 2
Ballycoolin; Dublin 15, Ireland
Phone: +353-18975900; Fax: +353-18293915
E-mail: medical@topcon.ie; www.topcon.ie



TOPCON CORPORATION

75-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8580, Japan.
Phone: 3-3558-2523/2522, Fax: 3-3960-4214, www.topcon.co.jp